



UTERINE
MANIPULATOR



VECTEC®
SINGLE USE
LAPAROSCOPIC
DEVICES



VECTEC®

UTERINE MANIPULATOR

Original design using suction
to get a complete uterine control




Reference	Designation	Diameter	Length	Packaging
CAUT60	Uterine manipulator	6 mm	90 mm	10 units per box

The uterine manipulator set consists in:

- A soft and atraumatic uterine dilator, double ended to perform hysteroscopy. The small size of the dilator preserves the cervix and the endometrium.
- The manipulator
 - an ergonomic handle design enables an optimal uterine manipulation during both diagnostic and surgical procedures
 - the atraumatic perforated blunt tip provides a perfect upholding of the cervix and the uterine cavity on the manipulator, thanks to the suction effect.
 - the stainless steel rod of the manipulator can be bended to obtain an adapted curved shape, if needed.
- Two reducers to be added or not according to the hysteroscopy.
- A suction tubing with Y shaped connector to be connected to the vacuum installation.

A fertility test with methylene blue can be carried out with the uterine manipulator.

The suction effect is obtained either through the vacuum installation and the aspiration tube or with a 50ml syringe.

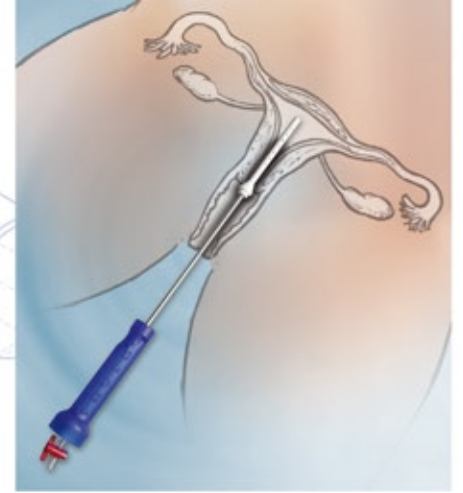

MD Class IIa
Sterilized by ETO

See instructions for use

CE
0459

G
MED
NF EN ISO 13485



International patent N° EP/2231040 and WO/2009/074844.

Manufactured by:

Péters Surgical
42, rue Benoit Frachon
93013 Bobigny Cedex - FRANCE
Tel: +33 (0)1 48 10 62 62
Fax: +33 (0)1 48 10 22 99

commercial@peters-surgical.com
www.peters-surgical.com

Péters
SURGICAL

98L71-GB / 05-2018

Non contractuel patient. Graphisme : Poissonniers - Octobre 2015